


H3K79me3 polyclonal antibody

Catalog # PAB14122 Size 50 ug


Applications

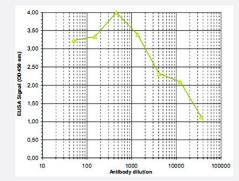

Enzyme-linked Immunoabsorbent Assay

ELISA was performed using a serial dilution of H3K79me3 polyclonal antibody (Cat # PAB14122) in antigen coated wells. The antigen used was a peptide containing the histone modification of interest.

By plotting the absorbance against the antibody dilution, the titer of the antibody was estimated to be 1 : 10,000.

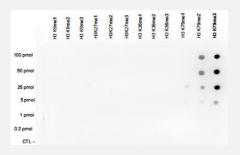

Dot Blot

To determine the cross reactivity of H3K79me3 polyclonal antibody (Cat # PAB14122) with other histone H3 modifications, a Dot Blot analysis was performed.

Tested histone H3 modifications include mono- and dimethylation of the same lysine and mono-, di- and trimethylation of lysine 9, 27 and 36. One hundred to 0.2 pmol of peptide containing the respective histone modification were spotted on a membrane and detected with the antibody diluted 1:50,000.

Figure shows a high specificity for the modification of interest.

Specification	
Product Description	Rabbit polyclonal antibody raised against synthetic peptide of H3K79me3.
Immunogen	A synthetic peptide (conjugated with KLH) corresponding to region of histone H3 containing the trime thylated lysine 79 (H3K79me3).
Host	Rabbit
Reactivity	Human, Mouse
Form	Liquid



Product Information

Recommend Usage	ELISA (1:500-1:1000) Dot Blot (1:50000) Western Blot (1:500) ChIP (5 ug/ChIP) The optimal working dilution should be determined by the end user.
Storage Buffer	In PBS (0.05% sodium azide, 0.05% proclin 300)
Storage Instruction	Store at -20°C. For long term storage store at -80°C. Aliquot to avoid repeated freezing and thawing.
Note	This product contains sodium azide: a POISONOUS AND HAZARDOUS SUBSTANCE which shoul d be handled by trained staff only.

Applications

- ChIP
- Western Blot
- Enzyme-linked Immunoabsorbent Assay

ELISA was performed using a serial dilution of H3K79me3 polyclonal antibody (Cat # PAB14122) in antigen coated wells. The antigen used was a peptide containing the histone modification of interest.

By plotting the absorbance against the antibody dilution, the titer of the antibody was estimated to be 1:10,000.

Dot Blot

To determine the cross reactivity of H3K79me3 polyclonal antibody (Cat # PAB14122) with other histone H3 modifications, a Dot Blot analysis was performed.

Tested histone H3 modifications include mono- and dimethylation of the same lysine and mono-, di- and trimethylation of lysine 9, 27 and 36. One hundred to 0.2 pmol of peptide containing the respective histone modification were spotted on a membrane and detected with the antibody diluted 1:50,000.

Figure shows a high specificity for the modification of interest.

Publication Reference

Dynamic regulation of histone lysine methylation by demethylases.

Shi Y, Whetstine JR.

Molecular Cell 2007 Jan; 25(1):1.





• Intra- and inter-nucleosomal protein-DNA interactions of the core histone tail domains in a model system.

Zheng C, Hayes JJ.

The Journal of Biological Chemistry 2003 Apr; 278(26):24217.

Translating the histone code.

Jenuwein T, Allis CD.

Science 2001 Aug; 293(5532):1074.